

Pin Dot



Dear Seating Professional:

When you order a seating system, you want three things: (1) a quality product, (2) reliable service, and (3) solid information. At Pin Dot we provide all three.


Our seating systems support clients in a functional and comfortable position. The custom-molded Contour-U cushions are guaranteed to fit even your toughest clients, and they are easily mounted on the Impala mobility base. For mildly involved wheelchair users EndoFlex is the best seat. It fits both manual and powered wheelchairs.

All products are backed by Pin Dot's unmatched service. If a cushion doesn't fit your client exactly as you want it to, we'll re-make it at no extra charge. We'll even take the cushion back and give you a full refund, if you feel it just isn't quite right.

Pin Dot provides all the information you need. We conduct three-day training courses, advanced training courses, and annual seating seminars. In addition, we offer a complete line of educational videotapes and publish two newsletters.

If you care about seating as much as we do, give us a call and let's work together.

Michael Silverman, C.O.
President

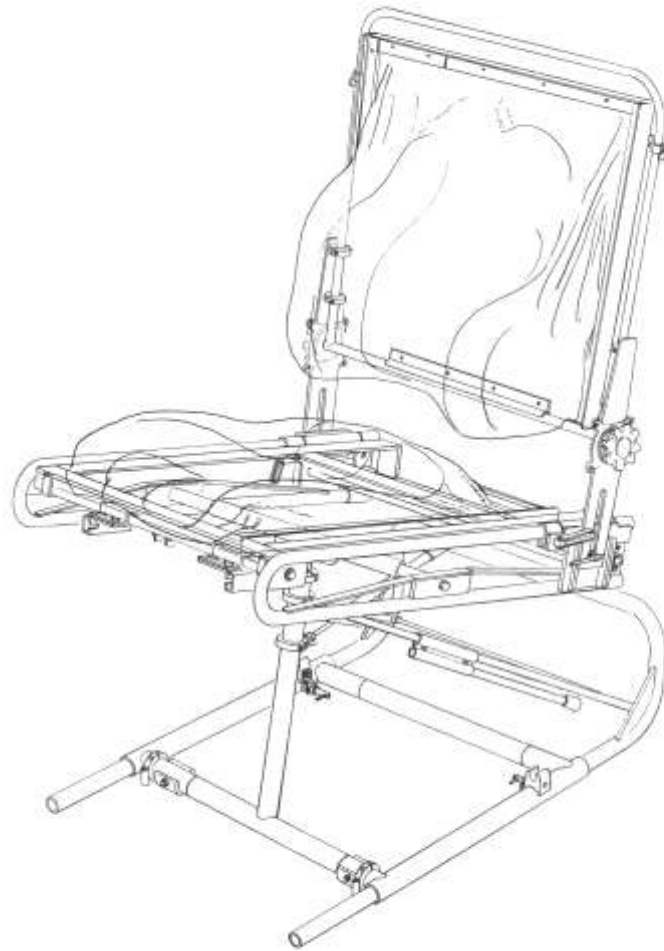

Simulator

1 Endoskeletal Design

The unique, endoskeletal (internal framework) design makes the hardware invisible. Nylon strapping provides suspension, flexibility, and support so EndoFlex never bottoms out.

2 Comfortable Contours

The seat contours support your hips and thighs comfortably. The anatomically-designed back follows your spine's curve and supports it in its natural position.

3 Dynamic Back

The dynamic lumbar support allows for free and easy movement. You can lean back, stretch, or shift weight. It's even easier to propel your wheelchair.

If you prefer a fixed back, the lumbar support can be locked in the position that is most comfortable for you.

4 Custom Fit

The seat is depth-adjustable, giving added support under your hips and thighs. The back is height-adjustable so the dynamic lumbar support fits precisely in your lower back.

The back-to-seat angle can be changed so your body is supported in the most comfortable and functional position.

The side supports are bendable and can be sloped to provide comfortable support.

5 Luxurious Cover

The seat and back are individually molded and covered with luxurious fabric. The fabric is soft, yet extremely durable. It will stand up to years of heavy use. You can even wash it with soap and water.

6 Easy Installation

Using patented mounting brackets, your dealer can install EndoFlex in just a few minutes. There is no need to drill holes in the frame or cut it, so your wheelchair remains completely intact.

7 Ergonomic Fit


Unlike old-fashioned cushions, EndoFlex does not raise you up or push you forward. The seat remains flush with the seat rails and ensures an ergonomic fit. The wheels are kept within easy reach, and your knees continue to fit under the table.

8 Portable Design

EndoFlex is easy to take with you. Just remove the Clevises, fold the seat in half, and lift it by the built-in handle.

9 Available Options

EndoFlex comes in many different styles and colors. Vinyl upholstery is also available. Contact your local dealer for more information.

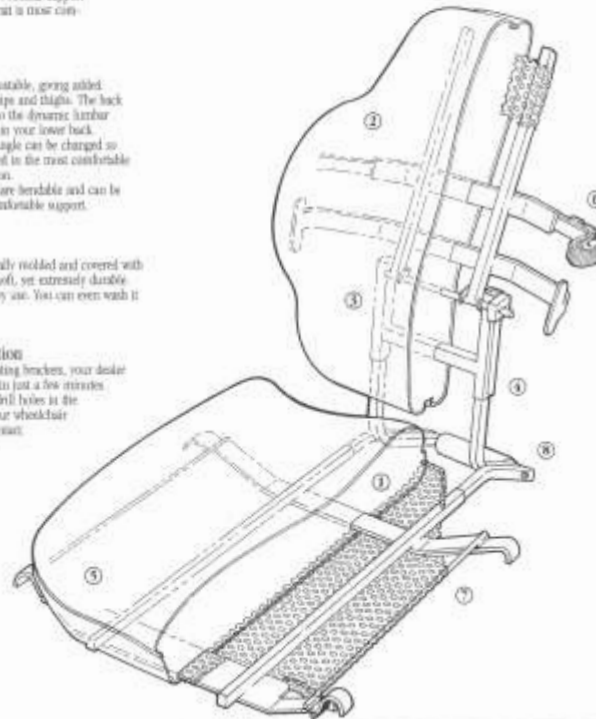

We've given rehab seating a whole new angle.

A mobility base that's easy to use is not news. But a mobility base that's easy to use, grows with the user, and *tilts*? That's news. And that's Impala.

Impala's dual cylinder tilt adjustment lets you adjust the angle of your seat plus or minus 30 degrees. That one small feature can help minimize the effects of hypotonia, fixed spinal deformities, or primitive flexor tone. And since the seat pivot is under the knees, you can tilt to the max and still fit under tables and fixtures.

If you value cost efficiency, you'll be pleased to discover it's hard to improve an Impala. Our base can handle up to 250 pounds, expand from 14 to 18 inches wide, and accommodate your choice of cushions, including Gortone-4.

Impala also lets you choose the wheel that's right for you: 12-inch, for getting in and out of tight spots on level surfaces; 24-inch, for smoothing over bumps and dips; or our optional 20-inch, which brings Impala's seat rail height to just 17 inches for superior foot/ground contact. And Impala's opens front legs let you position legs for comfort — even legs with tight harnesses — as easily as you can position them for mobility.

Is Impala hard to assemble? Not at all. Just install the wheels and you're ready to go: all the hardware you need is included. If you're building a seating system from the ground up, we've even put together an interfacing kit, complete with pins, dual links, armrests, foot rests, and push handles.



IMPALA™

The Impala mobility base.
It's a new slant on mobility.

Pin Dot Products
600 Gross Point Road
Niles, Illinois 60048-4027
Phone 800-451-5553



Impala (post Invacare acquisition)